

Архангельск (8182)63-90-72	Калининград (4012)72-03-81	Нижний Новгород (831)429-08-12	Смоленск (4812)29-41-54
Астана +7(7172)727-132	Калуга (4842)92-23-67	Новокузнецк (3843)20-46-81	Сочи (862)225-72-31
Белгород (4722)40-23-64	Кемерово (3842)65-04-62	Новосибирск (383)227-86-73	Ставрополь (8652)20-65-13
Брянск (4832)59-03-52	Киров (8332)68-02-04	Орел (4862)44-53-42	Тверь (4822)63-31-35
Владивосток (423)249-28-31	Краснодар (861)203-40-90	Оренбург (3532)37-68-04	Томск (3822)98-41-53
Волгоград (844)278-03-48	Красноярск (391)204-63-61	Пенза (8412)22-31-16	Тула (4872)74-02-29
Вологда (8172)26-41-59	Курск (4712)77-13-04	Пермь (342)205-81-47	Тюмень (3452)66-21-18
Воронеж (473)204-51-73	Липецк (4742)52-20-81	Ростов-на-Дону (863)308-18-15	Ульяновск (8422)24-23-59
Екатеринбург (343)384-55-89	Магнитогорск (3519)55-03-13	Рязань (4912)46-61-64	Уфа (347)229-48-12
Иваново (4932)77-34-06	Москва (495)268-04-70	Самара (846)206-03-16	Челябинск (351)202-03-61
Ижевск (3412)26-03-58	Мурманск (8152)59-64-93	Санкт-Петербург (812)309-46-40	Череповец (8202)49-02-64
Казань (843)206-01-48	Набережные Челны (8552)20-53-41	Саратов (845)249-38-78	Ярославль (4852)69-52-93

Единый адрес для всех регионов: <http://unispectr.nt-rt.ru/> || urc@nt-rt.ru

Центрифуги PowerSpin™ EX

Технические характеристики

PowerSpin™ EX Centrifuge

PowerSpin™ EX Centrifuge With 24 place horizontal rotor
this centrifuge is a great answer for a busy laboratory.


The PowerSpin™ EX has variable speed 500 – 3500 RPM,
to insure a clean and complete separation. Quiet and
dependable, with a maintenance free motor, with stainless
steel tube shields it can be used to process blood, urine,
fecal and other specimens.

The all digital front panel with 3 pre-programmed options

- 3200 RPM for 10 minutes for Blood
- 1800 RPM for 5 minutes for Urine
- 1300 RPM for 10 minutes for Fecal

Custom setting are easy and intuitive with all digital panel
and oversized display.

Model #C8924

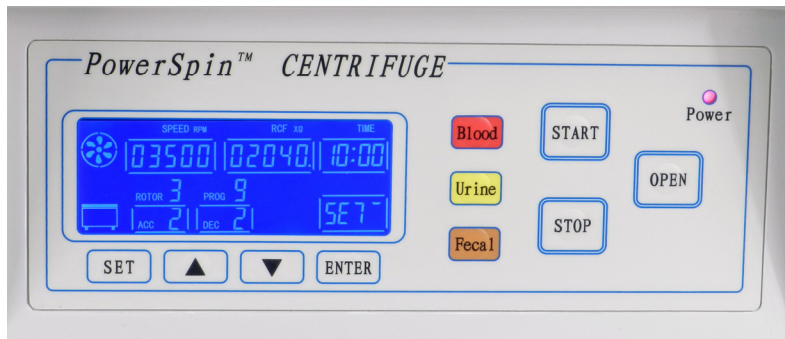

- Zero RPM Safety Latching System prevents lid from being opened while rotor is spinning
- Quiet running, only 50 decibels at 3500 RPM
- 99 minute digital timer, 1 sec. increment
- 500-3500 RPM, 100 RPM increment

- 3 pre-programmed options for blood, fecal and urine
- Remembers last setting for quick repeat runs
- Adjustable Acceleration and Deceleration
- Stainless steel rotor chamber for easy cleaning
- Available in **110V and 220V**

MFG Part #	Description	UOM
C8924	PowerSpin EX, 24 place horizontal spin centrifuge. 110V	EA
C8924E	PowerSpin EX, 24 place horizontal spin centrifuge. 220V	EA

Архангельск (8182)63-90-72	Калининград (4012)72-03-81	Нижний Новгород (831)429-08-12	Смоленск (4812)29-41-54
Астана +7(7172)727-132	Калуга (4842)92-23-67	Новокузнецк (3843)20-46-81	Сочи (862)225-72-31
Белгород (4722)40-23-64	Кемерово (3842)65-04-62	Новосибирск (383)227-86-73	Ставрополь (8652)20-65-13
Брянск (4832)59-03-52	Киров (8332)68-02-04	Орел (4862)44-53-42	Тверь (4822)63-31-35
Владивосток (423)249-28-31	Краснодар (861)203-40-90	Оренбург (3532)37-68-04	Томск (3822)98-41-53
Волгоград (844)278-03-48	Красноярск (391)204-63-61	Пенза (8412)22-31-16	Тула (4872)74-02-29
Вологда (8172)26-41-59	Курск (4712)77-13-04	Пермь (342)205-81-47	Тюмень (3452)66-21-18
Воронеж (473)204-51-73	Липецк (4742)52-20-81	Ростов-на-Дону (863)308-18-15	Ульяновск (8422)24-23-59
Екатеринбург (343)384-55-89	Магнитогорск (3519)55-03-13	Рязань (4912)46-61-64	Уфа (347)229-48-12
Иваново (4932)77-34-06	Москва (495)268-04-70	Самара (846)206-03-16	Челябинск (351)202-03-61
Ижевск (3412)26-03-58	Мурманск (8152)59-64-93	Санкт-Петербург (812)309-46-40	Череповец (8202)49-02-64
Казань (843)206-01-48	Набережные Челны (8552)20-53-41	Саратов (845)249-38-78	Ярославль (4852)69-52-93

Единый адрес для всех регионов: <http://unispectr.nt-rt.ru/> || urc@nt-rt.ru